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PL

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EVOREI
Boost 2
USER'S MANUAL

EN

Dear Customer, Thank you for trusting our brand and buying Evorei BOOST 2. The product you have chosen is a device designed for maintaining oral hygiene. The device is not a medical product and is designed for household use only. **Before using the device, please follow the user manual carefully.** Failure to comply with the rules presented in the instructions can cause injuries or device damage. The manufacturer is not liable for the damages resulting from not following the instructions.

01. SAFETY

- Children 8 years and older and people with reduced motor, mental or cognitive abilities can use the device only under guardian supervision.
- It is not recommended to use the device longer than 5 minutes or more often than every 2 hours.
- Do not use the device if you have sores or cuts inside your mouth.
- Remove all jewelry from your mouth before using the device.
- Do not toss or put any objects into the device when open.
- If your doctor recommended you to take antibiotics before doing any dental work, consult your dentist before using the device.
- The device is not a toy.
- Cleaning and maintenance should not be performed by children without adult supervision.
- Do not disassemble the device except for when it is supposed to be utilized.
- Do not fill the tank with water over 40°C.
- Check the device before every use. If it appears to be damaged, do not use it.
- Do not hit or drop the device. If you drop or hit the device, and it becomes damaged – stop using it immediately.
- Store the charging cord away from heated surfaces.
- Do not power the device without previously mounting a jet tip.
- Do not immerse the device in water or other liquids.
- Keep the tip of the water flosser out of reach of children – choking hazard!
- Mouthwash can damage the device. Do not pour mouthwash inside the water flosser tank.
- If you feel pain or discomfort during or after irrigation, or fear of experiencing other reactions to the water flosser, consult your doctor.
- Do not press the tip of the flosser too hard against your gums or teeth.
- Use only original tips included in the set.

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- The device is designed to clean the spaces between teeth and tongue and to improve gum health. Do not use the device for other reasons than what stated in the user manual!
- Do not direct the water stream directly under the tongue, inside the ear, nose, and other delicate places. The device can produce such pressure that can damage these sensitive spots.
- The device is not designed to be used by children and people with reduced mental or cognitive capabilities.
- Do not try and fix the device on your own, or modify and disassemble the device parts.
- Do not put the device in a place where it can fall into a tub or sink.


02. ORAL HYGIENE IRRIGATOR

Designed to clean the spaces between permanent braces, crowns, and dental bridges. It gets rid of up to 99% plaque and perfectly cleans the spaces between teeth and below gum lines.

Standard teeth cleaning: Washing your teeth with a traditional brush does not provide a thorough removal of all the leftover food particles and bacteria from between the teeth.

The range of dental floss is also limited - it can reach only the upper gum line.

Irrigator will precisely remove food particles residue, and pathogenic bacteria from the spaces between teeth.

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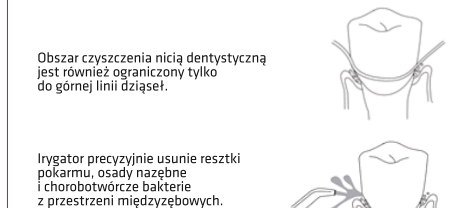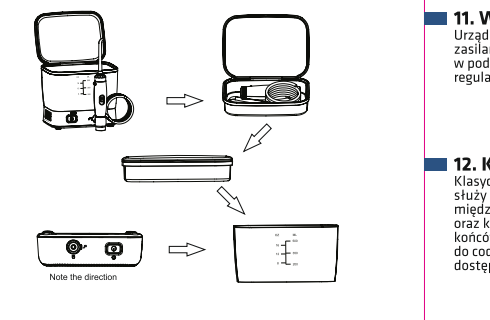
16. WARRANTY/ CLAIMS

This product is covered by a 24-month warranty. The warranty conditions are available on the website: <https://evorei.pl/warranty/>
Complaints should be submitted with the complaint form which can be found at <http://reklamacje.igktrade.pl/>
Contact details and the address of the service can be found on the website: www.evorei.pl

The crossed out trash can symbol indicates that unusable electrical or electronic devices, its accessories (such as power supplies, cords) or components (for example batteries, if included) cannot be disposed of alongside with household waste. In order to dispose of the devices or its components (for example, batteries) deliver the device to the collection point, where it will be accepted free of charge. Disposal is subject to the recast version of the WEEE Directive (2012/19 / EU) and the Directive on batteries and accumulators (2006/66 / EC). Proper disposal of the device prevents degradation of the natural environment. Information about the collection points of the facilities is issued by the competent local authorities. Incorrect disposal of waste is subject to penalties provided for by the law in force in the given area.

16. IMPORTER:

KGK Trade sp. z o.o. sp. k.
Ujastek 5B
31-752 Kraków
Poland

Made in P.R.C.